

По вопросам продаж и поддержки обращайтесь:

Алматы (7273)495-231
Ангарск (3955)60-70-56
Архангельск (8182)63-90-72
Астрахань (8512)99-46-04
Барнаул (3852)73-04-60
Белгород (4722)40-23-64
Благовещенск (4162)22-76-07
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Владикавказ (8672)28-90-48
Владимир (4922)49-43-18
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89

Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Иркутск (395)279-98-46
Казань (843)206-01-48
Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Коломна (4966)23-41-49
Кострома (4942)77-07-48
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Курган (3522)50-90-47
Липецк (4742)52-20-81

Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41
Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Ноябрьск (3496)41-32-12
Новосибирск (383)227-86-73
Омск (3812)21-46-40
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Петрозаводск (8142)55-98-37
Псков (8112)59-10-37
Пермь (342)205-81-47

Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Саранск (8342)22-96-24
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78
Севастополь (8692)22-31-93
Симферополь (3652)67-13-56
Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Сургут (3462)77-98-35
Сыктывкар (8212)25-95-17
Тамбов (4752)50-40-97
Тверь (4822)63-31-35

Тольятти (8482)63-91-07
Томск (3822)98-41-53
Тула (4872)33-79-87
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Улан-Удэ (3012)59-97-51
Уфа (347)229-48-12
Хабаровск (4212)92-98-04
Чебоксары (8352)28-53-07
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Чита (3022)38-34-83
Якутск (4112)23-90-97
Ярославль (4852)69-52-93

Россия +7(495)268-04-70

Казахстан +7(7172)727-132

Киргизия +996(312)96-26-47

eus@nt-rt.ru | <https://escolifesciences.nt-rt.ru>

Laminar Flow Storage Cabinet

Introduction

Esco laminar flow cabinets are the premium selection for the discerning user. It offers a combination of value, high quality construction, low operating noise levels, and a wide product range to suit all budgets, from an industry leader.

Basic Principles

Esco laminar flow storage cabinets make a positive contribution to maintaining the cleanliness of a cleanroom environment.

- ULPA-filtered airflow keeps garments clean during storage and handling.
- High quality ULPA filters which utilizes an improved mini-pleated separation technique to maximize surface area improving efficiency and extending the filter life. Filters operate at a typical efficiency of >99.999% at 0.1 to 0.3 micron sizes, providing superior product protection over conventional HEPA filters.
- Enables garments to be stored in a visible and organized manner.
- Esco laminar flow cabinets provide ISO Class 3 air cleanliness within the work zone as per ISO 14644.1, 100 times cleaner than the usual Class 5 classification.
- An additional disposable pre-filter on all models traps large particles in the inflow air prior to reaching the main filter, protecting it against damage and prolonging its life.

Standard Features


- Reliable rocker switches control the fan and lights and a Minihelic™ pressure gauge monitors cabinet operation.

- Built-in warm white, electronically ballasted, 5000k lighting provides excellent illumination of the work zone and reduces operator fatigue. The reliable lighting system is zero-flicker and instant start.
- All components are designed for maximum chemical resistance and enhanced durability for a long service life.
- The main body of the cabinet is constructed with industrial-grade electrogalvanized steel.
- The cabinet can be designed as mobile with caster wheels or static via built-in leveling feet.
- All cabinet components are clean room compatible. Isocide eliminates 99.9% of surface bacteria within 24 hours of exposure.

Enhanced Filtration System

- Esco laminar flow cabinets provide ISO Class 3 air cleanliness within the work zone as per ISO 14644.1, 100 times cleaner than the usual Class 5 classification on cabinets offered by the competition.
- High quality ULPA filters utilizing an improved mini-pleated separation technique to maximizes surface area improves efficiency and extends the filter life. Filters operate at a typical efficiency of >99.999% at 0.1 to 0.3 micron sizes, providing superior product protection over conventional HEPA filters.

Cabinet Filtration System

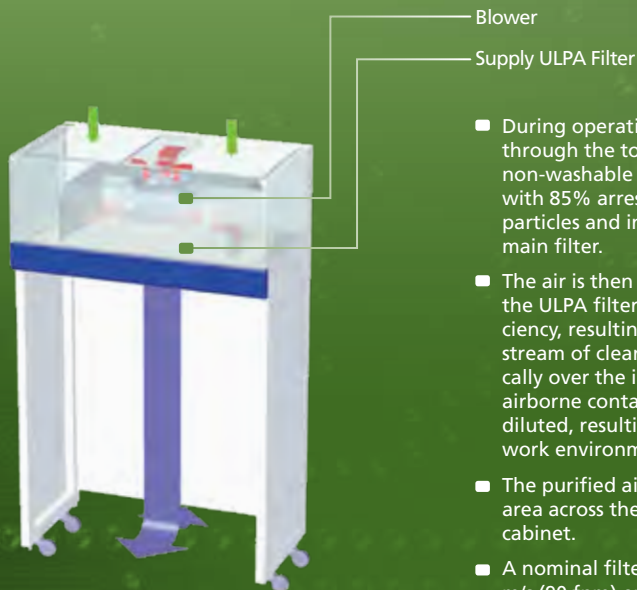

- During operation, room air is drawn through the top of the cabinet via a non-washable polyurethane pre-filter with 85% arrestance, trapping larger particles and increasing the life of the main filter.
- The air is then forced evenly through the ULPA filter with >99.999% efficiency, resulting in a unidirectional stream of clean air projected vertically over the internal work zone. All airborne contaminants are flushed and diluted, resulting in a particulate-free work environment.
- The purified air then leaves the storage area across the entire open front of the cabinet.
- A nominal filter face velocity of 0.45 m/s (90 fpm) ensures that there is a sufficient amount of air changes within the enclosed area of the cabinet in order to maintain cleanliness.
- Esco laminar flow storage cabinets incorporate permanently lubricated direct drive centrifugal blowers.
- The energy efficient external rotor motor design reduces operating costs, noise, and vibration levels.
- Built-in solid state variable speed controllers with integral noise filters offering flexible adjustment from zero to maximum setting.
- Each cabinet is individually factory tested for safety and performance in compliance with international standards.
- All electrical components are UL listed or UL recognized, ensuring superior electrical safety.
- All Esco laminar flow storage cabinets meet general safety requirements set by independent testing laboratories.

- ULPA-filtered air
- Unfiltered / potentially contaminated air
- Room air / Inflow air

2


	Cabinet Performance	Air Quality	Filtration	Electrical Safety
Standards Compliance	EN 12469 IEST-RP-CC002.2, Worldwide	ISO 14644.1, Class 3, Worldwide AS 1386 Class 1.5, Australia JIS B9920 Class 3, Japan	EN-1822 (H14), Europe IEST-RP-CC001.3, Worldwide IEST-RP-CC007.1, Worldwide IEST-RP-CC034.1, Worldwide	IEC 61010-1, Worldwide EN 61010-1, Europe UL 61010-1, USA CAN/CSA-22.2, No.61010-1



Mini-pleat Separatorless Filter (left) vs. Conventional Aluminium Separator Filter (right)



Esco cabinets use Swedish Camfil Farr® mini-pleat filters without aluminum separators to increase filter efficiency, minimize the chance of leakage, and to prolong filter life. Filters include a lightweight aluminum frame for structural stability and elimination of swelling common to conventional wood frames.


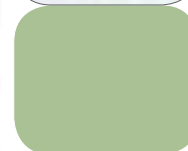


General Specifications, Laminar Flow Storage Cabinet

Model		EQU/04-LFSC	EQU/06 LFSC
Nominal Size		1.2 meters (4')	1.8 meters (6')
External Dimensions (W x D x H)		1340 x 740 x 2240 mm 52.8" x 29.1" x 88.2"	1950 X 740 X 2240 mm 76.8" x 29.1" x 88.2"
Internal Storage Area, Dimensions (W x D x H)		1240 x 645 x 1790 mm 48.8" x 25.4" x 70.5"	1850 x 645 x 1790 mm 72.8" x 25.4" x 70.5"
Storage Capacity		16 garments on hangers (4' model)	24 garments on hangers (6' model)
Average Airflow Velocity		0.45 m/s (90 fpm)	
Pre-Filter		Disposable and non-washable polyester fibers with 85% arrestance / EU3 rated	
ULPA Filter Typical Efficiency		99.999% for particles size at 0.3 microns	
Sound Emission Per IEST-RP-CC002.2		61 dBA	63 dBA
Fluorescent Lamp Intensity At Zero Ambient		>800 Lux (74 foot candles)	
Cabinet Construction	Main Body	1.2mm (0.05") 18 gauge electro-galvanised steel with white oven-baked epoxy powder-coated finish	
Electrical	220-240V, AC, 50Hz, 1Ø	EQU/04-LFSC	EQU/06-LFSC
	Cabinet Nominal Power	378 W	628 W
	Cabinet Full Load Amps (FLA)	1.8 A	4 A
	Cabinet BTU	1290	2143
Net Weight		150 kg (331 lbs)	220 kg (484 lbs)
Gross Weight		229.5 kg (506 lbs)	311.6 kg (687 lbs)
Shipping Dimensions, Maximum (W x D x H)		2150 x 950 x 1610 mm 84.6" x 37.4" x 63.4"	2150 x 950 x 1610 mm 84.6" x 37.4" x 63.4"
Shipping Volume, Maximum		3.29 m³ (116 cu.ft.)	3.29 m³ (116 cu.ft.)

3

Model LFSC Esco Laminar Flow Storage Cabinet

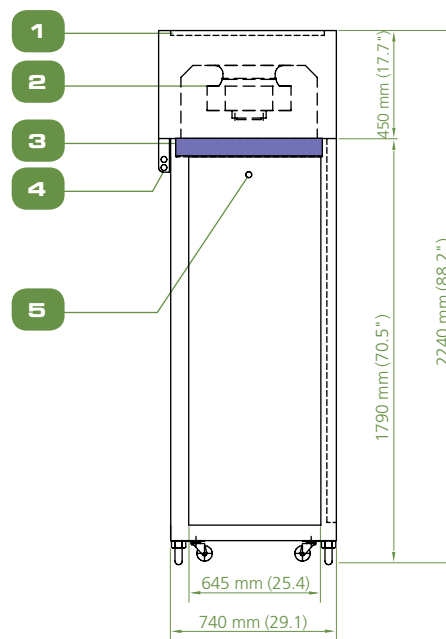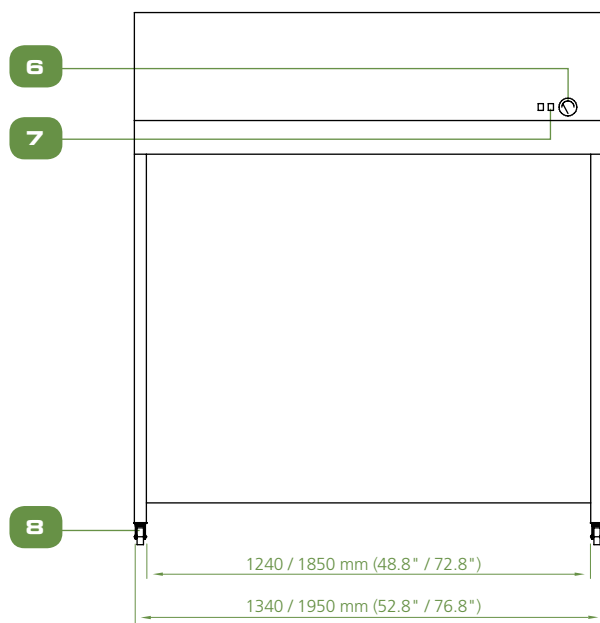

1. Pre-filter
2. Blower
3. ULPA filter
4. Fluorescent lamps
5. Stainless steel rod

6. Pressure gauge
7. Operating switches
8. Castors

По вопросам продаж и поддержки обращайтесь:

Алматы (7273)495-231
Ангарск (3955)60-70-56
Архангельск (8182)63-90-72
Астрахань (8512)99-46-04
Барнаул (3852)73-04-60
Белгород (4722)40-23-64
Благовещенск (4162)22-76-07
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Владикавказ (8672)28-90-48
Владимир (4922)49-43-18
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89
Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Иркутск (395)279-98-46
Казань (843)206-01-48

Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Коломна (4966)23-41-49
Кострома (4942)77-07-48
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Курган (3522)50-90-47
Липецк (4742)52-20-81
Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41
Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Ноябрьск (3496)41-32-12
Новосибирск (383)227-86-73

Омск (3812)21-46-40
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Петрозаводск (8142)55-98-37
Псков (8112)59-10-37
Пермь (342)205-81-47
Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Саранск (8342)22-96-24
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78
Севастополь (8692)22-31-93
Симферополь (3652)67-13-56
Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Сургут (3462)77-98-35

Сыктывкар (8212)25-95-17
Тамбов (4752)50-40-97
Тверь (4822)63-31-35
Тольятти (8482)63-91-07
Томск (3822)98-41-53
Тула (4872)33-79-87
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Улан-Удэ (3012)59-97-51
Уфа (347)229-48-12
Хабаровск (4212)92-98-04
Чебоксары (8352)28-53-07
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Чита (3022)38-34-83
Якутск (4112)23-90-97
Ярославль (4852)69-52-93

Россия +7(495)268-04-70

Казахстан +7(7172)727-132

Киргизия +996(312)96-26-47

eus@nt-rt.ru | | <https://escolifesciences.nt-rt.ru>